

# SunTech Vet BP Product Comparison

# SUNTECH

Utilizing an animal-specific blood pressure (BP) algorithm with our motion-tolerant technology, SunTech is dedicated to providing you with the most advanced solutions for reliably and accurately assessing animal BP and vital signs – an important advancement for pets and vets!

Our product line allows you to choose the BP and vital signs solution that works for your veterinary practice. Whether you need BP for wellness check-ups or a multiparameter monitor during procedures, we can find the solution for you.

The Vet20 is the simplest way to get a BP on an awake cat or dog. The Vet25 is designed to be flexible– for measuring BP in either the exam room or the operating room. The Vet30 combines our trusted BP technology with oxygen saturation and temperature monitoring for animals pre-, post-, or during a procedure.



**Vet20**



**Vet25**




**Vet30**

## Features

|   |                  |                     |                     |
|---|------------------|---------------------|---------------------|
| SunTech Advantage™ Motion Tolerant BP         | ✓                | ✓                   | ✓                   |
| Selectable BP Averaging                       | ✓                | ✓                   | ✓                   |
| Protective Armor Available in 3 Bright Colors | ✓                | ✓                   | ✓                   |
| Editable Clinical Alarms                      | ✓                | ✓                   | ✓                   |
| Universal AC Adapter                          | ✓                | ✓                   | ✓                   |
| Companion Animal Modes                        | ✓                | ✓                   | ✓                   |
| Equine Mode Options                           |                  | ✓                   | ✓                   |
| Bluetooth Capability for Data Transfer        |                  | ✓                   | ✓                   |
| PC Application for Reports and Data Trending  |                  | ✓                   | ✓                   |
| Interval BP Mode                              |                  | ✓                   | ✓                   |
| STAT BP Mode                                  |                  | ✓                   | ✓                   |
| Accuvet™ Temperature Measurement              |                  |                     | ✓                   |
| Accuvet™ or Masimo™ SpO2 Measurement          |                  |                     | ✓                   |
| Focus Screens                                 |                  |                     | ✓                   |
| Memory Storage                                | 50               | 960                 | 960                 |
| Battery Type                                  | AA Alkaline      | Rechargeable Li-Ion | Rechargeable Li-Ion |
| Weight  | 1.3 lb (0.59 kg) | 2.1 lb (0.95 kg)    | 2.65 lb (1.2 kg)    |
| Sizes of Veterinary Cuffs                     | #1 - #7          | #1 - #9             | #1 - #9             |

# SunTech Vet BP Products

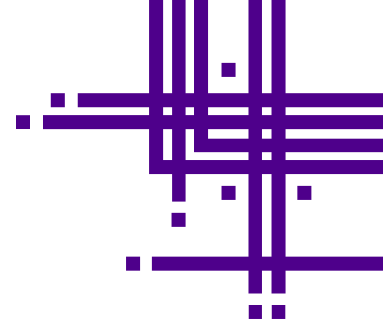

## Vet20

### Spot-Check Blood Pressure Monitor for Companion Animals

The portable SunTech® Vet20™ quickly and easily provides reliable blood pressure (BP) measurements on awake companion animals, ensuring a simpler assessment process for veterinarians and their staff members. Utilizing SunTech's trusted Advantage™ VET BP technology, the Vet20 quietly performs motion-tolerant, oscillometric BP measurements, delivering accurate results. The intuitive touch screen interface includes all of the necessary features for reviewing and analyzing BP data. Alarms can be silenced to avoid frightening the animal and saved measurement data is always retained when the device is powered off.

## Vet25 & Vet25E

### Spot-Check and Continuous BP Monitor for Animals

The portable and rechargeable SunTech® Vet25™ monitor quickly and easily provides reliable blood pressure (BP) measurements for wellness exams, spot-checks, and while monitoring companion animals pre- and post-procedure. Intervals can be set to automatically take BP readings during a monitoring period. Utilizing SunTech's trusted Advantage™ VET BP technology, the Vet25 quietly performs motion-tolerant BP measurements, delivering accurate results. The intuitive, touchscreen interface provides users all of the necessary features for reviewing and analyzing BP data. The Vet25 is also available as the Vet25E, which comes with an equine-specific blood pressure mode in addition to companion animal modes.



## Vet30 & Vet30E

### Continuous Vital Signs Monitor for Animals

The portable and rechargeable SunTech® Vet30™ monitor quickly and easily provides reliable blood pressure (BP), oxygen saturation (SpO2), and temperature measurements while monitoring companion animals pre- and post-procedure. Intervals can be set to automatically take BP readings during a monitoring period. Utilizing SunTech's trusted Advantage™ VET BP technology, the SunTech Vet30 quietly performs motion-tolerant BP measurements, delivering accurate results. In addition, the device provides continuous measurement of temperature and SpO2. The intuitive, touch screen interface provides users all the necessary features for reviewing, analyzing, and transferring vital signs data. The Vet30 is also available as the Vet30E, which comes with an equine-specific blood pressure mode in addition to companion animal modes.

## SunTech Veterinary Cuffs

Our veterinary blood pressure cuff was specifically created with pets and vets in mind and is compatible with any veterinary blood pressure monitor. From your tiniest ball of fur to your largest Clydesdale, our veterinary blood pressure cuffs are designed and sized to ensure accurate measurement for all of your 4-legged patients. The cuffs are made from soft, latex-free materials and feature rounded corners to provide a gentle cuff solution. Our vet BP cuffs are color-coded for quick size selection and include range indicators to ensure accurate cuff sizing and fit. SunTech Vet BP cuffs can be used with other veterinary BP monitors! Ask about our SunTech Vet Cuff Conversion Kit.

