

ABPM



SunTech[®] *Oscar 2*[™]

24-HR ABPM



Available with:



24-Hour Ambulatory Blood Pressure Monitoring System

SunTech Oscar 2

The Oscar 2™ system from SunTech Medical® sets the standard in 24-hour ambulatory blood pressure monitoring (ABPM) by utilizing proprietary motion-tolerant technology, meeting the highest clinical standards, optimizing patient comfort through innovation, giving clinicians practical interpretive analysis and incorporating the latest monitoring advances.

The Oscar 2 system is independently tested to meet the accuracy and performance requirements of the British Hypertension Society, the European Society of Hypertension International and the AAMI-SP10 standards. This system comes standard with both the patented Orbit™ ABPM cuff and interpretive software AccuWin™ Pro 4 and NOW features options for central blood pressure measurement and Bluetooth-connected True24 patient diary mobile app.

Standard Features Configuration Options

	ABPM	SphygmoCor	Bluetooth
SunTech Motion Tolerance	✓		
Clinical Validations	✓		✓
Patented Orbit ABPM Cuffs	✓	✓	
AccuWin Pro 4 software	✓	✓	✓

Orbit™ Blood Pressure Cuffs


SunTech's patented Orbit cuff is specifically engineered for ambulatory blood pressure monitoring applications. This cuff is designed with a soft sleeve that gently hugs the arm to maintain cuff placement and provide the ultimate in patient comfort for an extended period of time.



AccuWin Pro™ 4 Software

AccuWin Pro 4 is a user-friendly, Windows® based program that allows for configuration, analysis, interpretation and reporting of ABPM studies. You have the ability to rapidly decipher data, creating interpretive summaries, according to AHA, ESH, or JNC 7 guidelines, for:


- automated analysis of ABP levels for normality
- white coat syndrome
- overnight dipper status

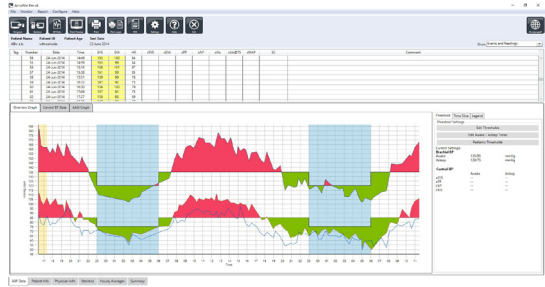

Optional Features

SphygmoCor Inside: SunTech Medical has partnered with AtCor Medical, the Gold Standard in non-invasive measurement of central BP and arterial stiffness, to bring you the first and only ABPM Monitor with SphygmoCor inside. The Oscar 2 ABPM system delivers a patient-focused solution with unmet comfort, data reliability and confidence.

Bluetooth with True24: The True24 patient diary mobile app* is designed to pair exclusively with the Oscar 2 ABPM device with optional Bluetooth. This easy-to-use app provides a tutorial on ambulatory blood pressure monitoring and prompts the patient to enter information about their activity, posture or symptoms experienced after a BP measurement is taken. This diary information is extremely valuable, further assisting clinicians with interpreting the ambulatory blood pressure study data and making hypertension treatment decisions.

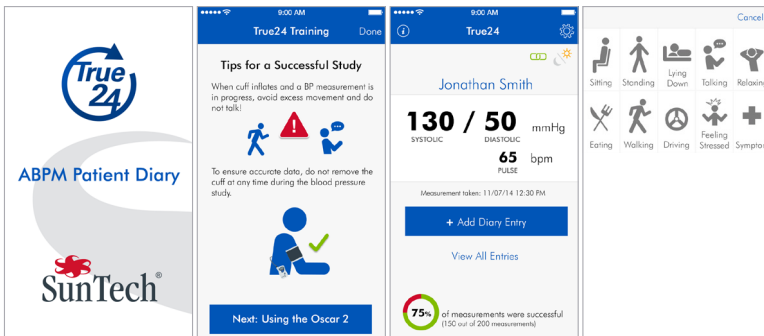

*True24 available only in the Apple iPhone App store

Technical Specifications

The Oscar 2 System Includes

- Oscar 2 monitor
- AccuWin Pro™ 4 (AWP4) software (Download)
- 2 Orbit™ ABPM cuffs (Adult and Adult Plus) Carrying pouch & patient belt
- USB interface cable
- User guide
- Standard 2 year warranty

True24 App for Bluetooth Systems


Available in the Apple App Store

General Specifications

- Validations:** BHS (A/A), ESH, ANSI/AAMI SP10:2002
- BP Technique:** Oscillometry with step deflation
- BP Range:** 25-260 mmHg
- HR Range:** 40-200 bpm
- Sampling Periods:** 3 adjustable (awake, asleep, and special)
- Time Intervals:** 5, 10, 15, 20, 30, 45, 60, 90 and 120 min.
- Event Key:** Program to allow additional patient readings
- Data Storage:** Flash memory stores up to 250 readings
- Power Source:** Two AA batteries
- Weight:** 284 grams, including batteries
- Dimensions:** 12 x 7 x 3 cm
- PC Requirements:** USB port, Windows® XP or later
- App Requirements:** Apple iPhone



With more than 30 years dedicated to developing our world-class blood pressure technology, we are confident that we don't miss a beat - so you won't either. It's our legacy. It's SunTech Medical. **It's The Difference In Clinical Grade.**


SunTech Medical, Inc.
507 Airport Boulevard, Suite 117
Morrisville, NC 27560-8200 USA
Tel: + 1.919.654.2300
1.800.421.8626
Fax: +1.919.654.2301
SunTechMed.com


EMERGO Europe
Prinsessegracht 20
2514 AP The Hague
The Netherlands

SunTech Medical (Shenzhen) Co., Ltd.
105 HuanGuan South Road, Suite 15 2~3/F
DaHe Community Guanlan,
LongHua District, Shenzhen
GuangDong PRC 518110
Tel.: + 86.755.29588810
+ 86.755.29588986 (Sales)
+ 86.755.29588665 (Service)
Fax: + 86.755.29588829
SunTechMed.com.CN





SunTech Oscar2 Part Numbers

Oscar 2 System Monitors

99-0133-00	Oscar 2 Monitor - Standard
99-0133-01	Oscar 2 Monitor - Bluetooth
99-0133-02	Oscar 2 Monitor - w/ SphygmoCor Inside
99-0133-03	Oscar 2 Monitor - w/ SphygmoCor Inside - Bluetooth
99-0133-10	Oscar 2 Monitor - Standard - Without AWP4
99-0133-11	Oscar 2 Monitor - Bluetooth - Without AWP4
99-0133-12	Oscar 2 Monitor - w/ SphygmoCor Inside - Without AWP4
99-0133-13	Oscar 2 Monitor - w/ SphygmoCor Inside - BT - Without AWP4
98-0056-00	Oscar 2™ PowerPack Accessory Kit

Orbit™ Comfort Blood Pressure Cuffs

98-0239-01	#1 Orbit™ ABPM Cuff (18-26 cm)
98-0239-02	#2 Orbit™ ABPM Cuff (26–34 cm)
98-0239-03	#3 Orbit™ ABPM Cuff (32–44 cm)
98-0239-04	#4 Orbit™ ABPM Cuff (42–55 cm)

Straps and Belts

98-0032-07	Oscar 2™ Pouch
98-0037-00	Oscar 2™ Belt

Extended Warranty Program

99-0027-00	One Year Extended Warranty
99-0027-03	Three Years Extended Warranty

Extended Warranty Terms and Conditions

The Extended Warranty must be purchased either:

- At the time of the original sale
- Prior to the end of the initial warranty period
- Prior to the end of the purchased extended warranty period

CPT Codes for ABPM Reimbursement

CPT code: 93784: Ambulatory blood pressure monitoring utilizing a system such as magnetic tape and/or computer disk, for 24 hours or longer; including recording, scanning analysis, interpretation, and report.

ICD-10 code: R03.0: White-coat hypertension. Elevated blood pressure reading without diagnosis of hypertension.