

24-Hour Ambulatory
Blood Pressure Monitoring System


SunTech[®] *Oscar2*[™]
24-HR ABP

Accurate BP Readings Day & Night



- + Highest Standards of Clinical Accuracy
- + Interpretive AccuWin Pro Software

- + Simple Set-Up & Operation
- + User-Centered Design

THE OSCAR 2™ SYSTEM FOR AMBULATORY BP MONITORING

SunTech Medical®, a company with a twenty-year history in blood pressure technology, presents the Oscar 2™ system for ambulatory blood pressure monitoring (ABPM). ABPM is considered the gold standard for measuring blood pressure by clinicians worldwide. It provides valuable diagnostic information that in-clinic and home blood pressure monitoring systems are incapable of measuring, including: bp variability, night-time bp, overnight dipping, and morning surge of bp.



+ Trusted Accuracy //

The Oscar 2 is the only ABPM available that is clinically validated to all three internationally recognized standards:

- British Hypertension Society (BHS) A/A
- European Society of Hypertension International Protocol (ESH)
- Association for the Advancement of Medical Instrumentation (AAMI) SP10


+ Designed for Patient Comfort & Compliance //

- Light and compact design
- Motion-tolerant technology reduces re-inflates and failed readings
- Dynamic inflation above previous systolic reduces measurement time
- Quiet pump operation
- Flash memory eliminates the possibility of lost data



+ Indications for Use^{1,2} //


- White-coat hypertension (reimbursable procedure)
- Drug resistant hypertension
- Masked hypertension
- Hypotensive symptoms with anti-hypertensive meds

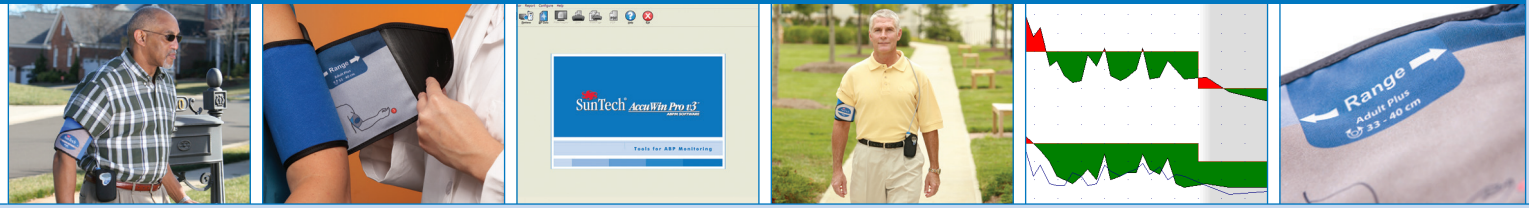

+ AccuWin Pro™ v3 //

AccuWin Pro™ v3 is a user-friendly, Windows® based program that allows maximum flexibility for configuration, analysis, interpretation and reporting of ABPM studies.

- **Interpretive Summary //** Automated analysis of ABPM data according to AHA, ESH or JNC 7 guidelines.
- **White-Coat Analysis //** The first hour of collected data is compared with the full 24-hour data to automatically indicate if the patient exhibits a white coat effect. Medicare reimbursement documentation is clearly displayed within the report.
- **Pediatric Threshold Calculator //** Pre-programmed pediatric ABP thresholds based on published guidelines improve workflow and analysis.
- **PDF Report Allows EHR/EMR Compatibility //** Create an un-editable, electronic report for storage in most EHR/EMR systems.
- **Single Page Report Option //** A simple one page report includes an ABP graph, interpretive summary, asleep dip analysis and blood pressure statistics with selected patient information.
- **Time-Slice Periods //** Provides statistical analysis on user-defined time windows.

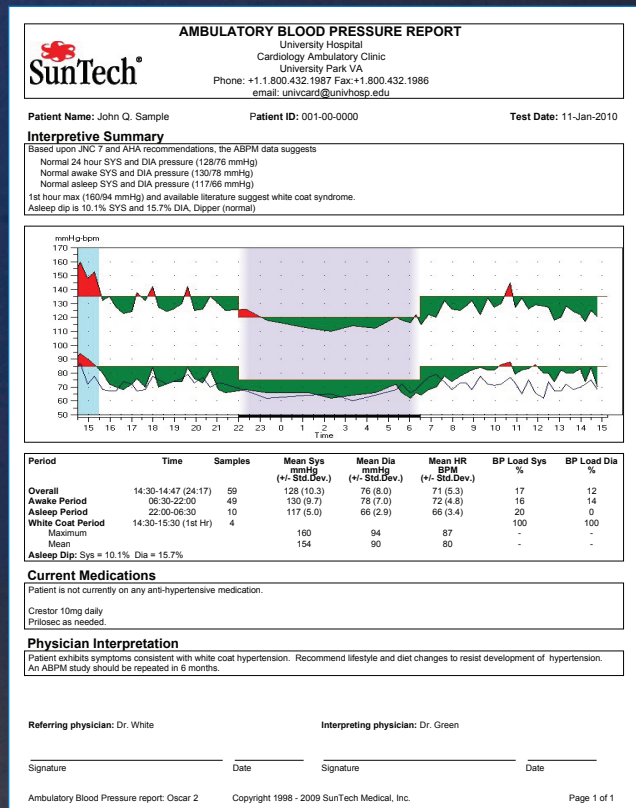

+ Orbit™ Blood Pressure Cuffs //

- Patented, stretch sleeve keeps the cuff in place both day and night
- Latex-free material
- Curved design fits larger arm sizes
- Removable bladder for machine washing





+ The Oscar 2™ System Includes:

- Oscar 2 ABP monitor
- AccuWin Pro™ v3 software (CD)
- 2 Orbit™ ABP cuffs (Adult and Adult Plus)
- Carrying pouch & patient belt
- USB interface cable
- User's guide
- Standard 2 year warranty

+ Specifications:

- Validations:** BHS (A/A), ESH, ANSI/AAMI SP10:2002
- BP Technique:** Oscillometry with step deflation
- BP Range:** 25-260 mmHg
- HR Range:** 40-200 bpm
- Sampling Periods:** 3 adjustable (awake, asleep, and special)
- Time Intervals:** 5, 10, 15, 20, 30, 45, 60, 90 and 120 min.
- Event Key:** Program to allow additional patient readings
- Data Storage:** Flash memory stores up to 250 readings
- Power Source:** Two AA batteries
- Weight:** 284 grams, including batteries
- Dimensions:** 12 x 7 x 3 cm
- PC Requirements:** Serial or USB port, Windows® XP or later



+ CPT Codes ABPM Reimbursement:

- | CPT 93784 | Description: |
|------------------|--|
| | ABPM, utilizing a system such as magnetic tape and/or computer disk, for 24 hours or longer, including recording, scanning analysis, interpretation, and report. |
| CPT 93786 | Recording only |
| CPT 93790 | Interpretation only |
| CPT 93788 | Scanning analysis with report |

Medicare covers ICD-9 code 796.2 which is described as elevated blood pressure reading without diagnosis of hypertension. Private insurance carriers may cover additional ICD-9 codes.

+ Product Numbers	+ Description
99-0012-00	Oscar 2 ABPM System
27-0037-A1	AccuWin Pro v3 software (CD)
98-0056-00	PowerPack Accessory Kit -
<small>* and highlighted items included in the Accessory Kit.</small>	Complete your system with a full line of Orbit cuffs and a 1-year extended warranty.
97-0090-02	PC USB Cable
98-0063-11 *	Small Adult Orbit Cuff (18-27 cm)
98-0063-12	Adult Orbit Cuff (25-35 cm)
98-0063-13	Adult Plus Orbit Cuff (33-40 cm)
98-0063-14*	Large Adult Orbit Cuff (39-46 cm)

+ Product Numbers	+ Description
98-0032-00	Oscar 2 Pouch
98-0037-00	Oscar 2 Belt
98-0036-00	Oscar 2 Shoulder Strap
82-0019-00	Patient Diaries (qty.50)
17-0007-00*	Battery Recharger + 4 NiMH AA batteries
80-0027-00	Oscar 2 English Manual
80-0026-00	AccuWin Pro v3 English Manual



SunTech Medical, Inc.
507 Airport Boulevard, Suite 117
Morrisville, NC 27560-8200 USA
Tel: + 1.919.654.2300
1.800.421.8626
Fax: + 1.919.654.2301

SunTech Medical, Ltd.
Oakfield Estate
Eynsham, Oxfordshire
England
Tel: +44 (0) 1865.884.234
Fax: +44 (0) 1865.884.235

SunTech Medical (Shenzhen) Co., Ltd.
105 HuanGuan South Road, Suite 15 2~3/F
DaHe Community Guanlan,
LongHua District, Shenzhen
GuangDong PRC 518110
Tel.: + 86-755-2958 8810
+ 86-755-29588986 (Sales)
+ 86-755-29588665 (Service)
Fax: + 86-755-2958 8829
www.SunTechMed.com.CN

1. The Seventh Report of Joint National Committee on Prevention, Detection, and Evaluation and Treatment of High Blood Pressure
2. World Health Organization - International Society of Hypertension