



 SunTech® **Vet25™**
BP MONITOR



Spot-Check and Continuous BP Monitor
for Companion Animals

SunTech Vet25

The portable and rechargeable SunTech® Vet25™ monitor quickly and easily provides reliable blood pressure (BP) measurements for wellness exams, spot-checks, and while monitoring companion animals pre- and post-procedure. Intervals can be set to automatically take BP readings during a monitoring period. Utilizing SunTech's trusted Advantage™ VET BP technology, the Vet25 quietly performs motion-tolerant BP measurements, delivering accurate results. The intuitive, touchscreen interface provides users all of the necessary features for reviewing and analyzing BP data.

Standard Features

Reliable: Due to SunTech's motion tolerant technology, the Vet25 has a high success rate for BP readings.

Simple to Use: Intuitive touch screen interface provides all necessary features for reviewing, analyzing, and transferring BP data. Designed for animal hospitals, not researchers.

Accurate: BP algorithm was specifically developed for companion animals by experts at a leading vet school.


Continuous Monitoring: Monitoring during procedures is simple with user-defined intervals for BP measurement.

Memory: Up to 960 data points are retained, even when device is powered off. Users can select specific measurements from the device memory and perform automated BP averaging.

Portable: Battery-operated device easily moved to different monitoring locations.

Rechargeable: Battery recharges when device is plugged into AC power, reducing disposable battery waste.

Data Transfer: Data can be transferred to a data trending and reporting software application for a PC using Bluetooth. Use the PC application to export reports in either PDF or CSV format.

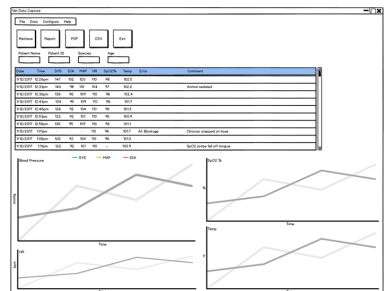



Technical Specifications

The Vet25 Includes:

- Vet25 System Monitor
- Protective Armour in Flamingo Pink, Peacock Blue, or Tree Frog Green
- Veterinary Blood Pressure Cuffs Sizes 1-6
- Patient Hose, 6 ft (1.8 m)
- AC Adapter

General Specifications

Systolic:	40 to 265 mmHg
Diastolic:	20 to 200 mmHg
Mean Arterial Pressure:	27 to 222 mmHg
Heart Rate:	25 to 300 BPM
Device Weight:	2.1 lbs (with battery), 953 grams (with battery)
Dimensions:	6.25" x 5" x 5.25", 15.9 cm x 12.7 cm x 13.3 cm
Warranty:	2 years on monitor 90 days on cuffs/hose

Accessories

SunTech Vet25 Armour

The removable protective armour is available in 3 colors and features a convenient handle.



Flamingo Pink



Peacock Blue



Tree Frog Green

AC Adapter

The Vet25 is a portable device that can easily travel from room to room when needed. Simply plug in the AC Adapter to recharge the lithium-ion battery in the Vet25. The AC Adapter includes 4 plug configurations to meet most international power requirements.

SunTech Veterinary Cuffs



Soft materials and rounded corners provide a gentle cuff solution. SunTech Medical BP cuffs are color-coded for quick size selection and range indicators ensure accurate cuff sizing and fit. Included with the Vet25 BP device are cuffs size 1–6.

SunTech Vet BP cuffs can be used with other veterinary BP monitors! Ask about our SunTech Vet Cuff Conversion Kit.



With more than 30 years dedicated to developing our world-class blood pressure technology, we are confident that we don't miss a beat - so you won't either. It's our legacy. It's SunTech Medical. ***It's The Difference In Clinical Grade.***



SunTech Medical, Inc.
507 Airport Boulevard, Suite 117
Morrisville, NC 27560-8200 USA
Tel: + 1.919.654.2300
1.800.421.8626
Fax: +1.919.654.2301
SunTechMed.com

SunTech Medical (Shenzhen) Co., Ltd.
105 HuanGuan South Road, Suite 15 2~3/F
DaHe Community Guanlan,
LongHua District, Shenzhen
GuangDong PRC 518110
Tel.: + 86.755.29588810
+ 86.755.29588986 (Sales)
+ 86.755.29588665 (Service)
Fax: + 86.755.29588829
SunTechMed.com.CN





SunTech Vet25 Part Numbers

Vet25 System Monitors

- 99-0170-00 Vet25 Monitor, Flamingo Pink Armour
- 99-0170-01 Vet25 Monitor, Peacock Blue Armour
- 99-0170-02 Vet25 Monitor, Tree Frog Green Armour

SunTech Veterinary BP Cuffs

- | | | | |
|----------------|----------------------------|-----------|----------------|
| 98-0400-80-VET | Cuff #1 | Box of 20 | Range 3-6 cm |
| 98-0400-81-VET | Cuff #2 | Box of 20 | Range 4-8 cm |
| 98-0400-82-VET | Cuff #3 | Box of 20 | Range 6-11 cm |
| 98-0400-83-VET | Cuff #4 | Box of 20 | Range 7-13 cm |
| 98-0400-84-VET | Cuff #5 | Box of 20 | Range 8-15 cm |
| 98-0400-F1 | Cuff #6 | Box of 20 | Range 12-19 cm |
| 98-0400-F3 | Cuff #7 | Box of 20 | Range 17-25 cm |
| 98-0240-00 | Cuffs 1-6 | Pack of 6 | |
| 91-0028-75 | Patient Hose, 6 ft (1.8 m) | | |
| 91-0028-76 | Patient Hose, 10 ft (3 m) | | |

Accessories

- 39-0195-03 Flamingo Pink Armour
- 39-0195-04 Peacock Blue Armour
- 39-0195-02 Tree Frog Green Armour
- 19-0020-00 AC Adapter
- 17-0026-00 Lithium-ion Battery Replacement

Extended Warranty Program

- 83-0022-00 One Year Extended Warranty
- 83-0022-01 Three Years Extended Warranty

Extended Warranty Terms and Conditions

The Extended Warranty can be purchased either:

- At the time of the original sale
- Prior to the end of the initial warranty period